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## WELCOME TO THE FAMILY

Congratulations on your purchase of the Candeon™ Rechargeable Oral Irrigator! This product combines the water pressure and pulsation generated by transmission pumps & pistons, which are used extensively in the automobile industry.

With high water pressure, the food debris lodged deeply between teeth and areas where brushing is hard to reach can be completely removed. With the pulsing water pressure, the irrigator gently massages your teeth and gums.

Thank you again for choosing Candeon™. Enjoy the experience of healthier teeth and a brighter smile!

This is a China patent product. All rights reserved by Candeon Technologies Ltd.

**CAUTION:**


**When using electrical appliances, basic safety precautions should always be followed to reduce the risk of fire, electric shock and/or injury to persons, including the following:**

- Read all the instructions before using the appliance.
- To reduce the risk of injury, close supervision is necessary when an appliance is used near children.
- Only use attachments recommended or sold by Candeon Technologies Limited or the authorized distributors.
- To reduce the risk of electric shock, do not put Candeon™ Rechargeable Oral Irrigator in water or other liquid. Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not handle charger with wet hands.
- Do not use the charger with a damaged charger cord or plug, or has been damaged in any manner.
- Do not place on or near a hot gas or electric burner or in a heated oven.
- Only use clean water or other solutions recommended by the dental professionals. Do not use mouthwashes/rinses that contain IODINE or Bleach.
- Do not use this appliance for other than intended use. Do not direct water under the tongue, into the ear, nose or other delicate areas that may cause injuries.
- Please read all directions in this manual carefully prior to operating this product.

**SAVE THESE INSTRUCTIONS.**

## PRODUCT DESCRIPTION

Product may vary slightly from what is illustrated.

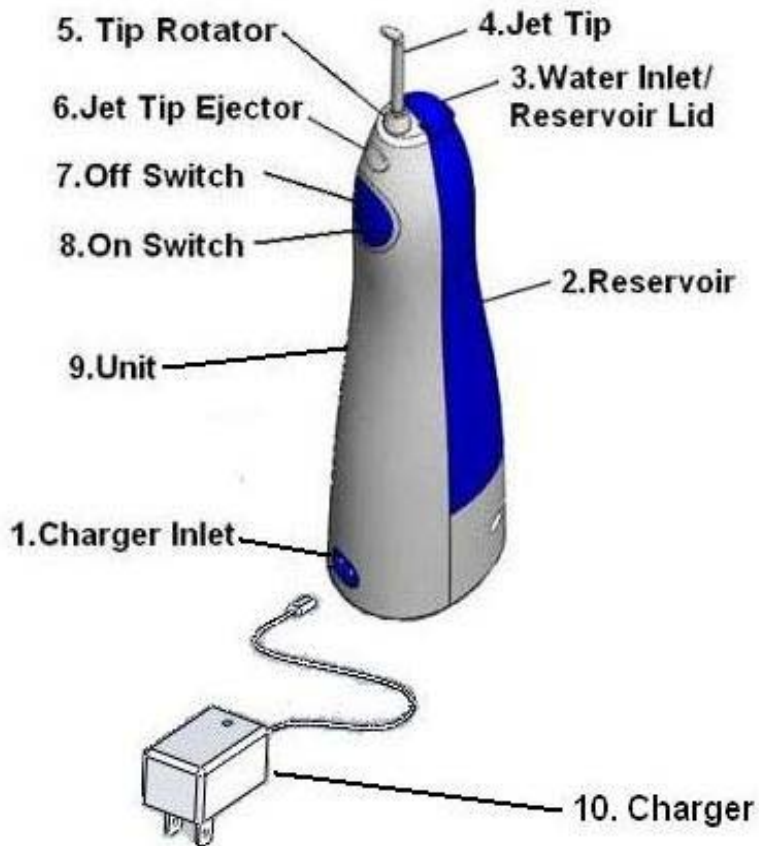1. Charger Inlet
2. Reservoir
3. Water Inlet/Reservoir Lid
4. Jet Tip
5. Tip Rotator
6. Jet Tip Ejector
7. Off Switch
8. On Switch
9. Unit
10. Charger


## GETTING STARTED

1. Read the Instruction Manual carefully.
2. Charge the Unit for 16 hours before the first use. **Important: The unit must be switched off before charging (A) !** The power LED indication on the Charger is lit when the power to the unit is on **(B)**.


**( A )**

3. Open the Reservoir Lid and fill the water into the Reservoir **( C )**.

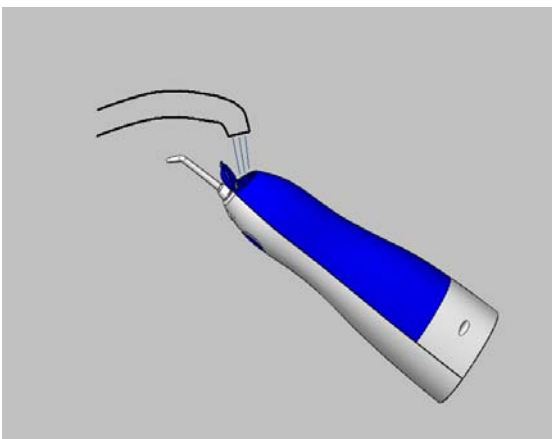
**( C )**

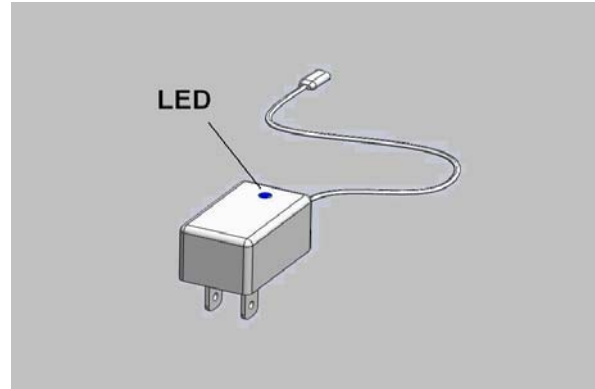**( B )**

4. Insert the Jet Tip into the center of the knob **( D )**. Note: There are 3 Tips with different colors in the original package which are prepared for the individual use of your family members.

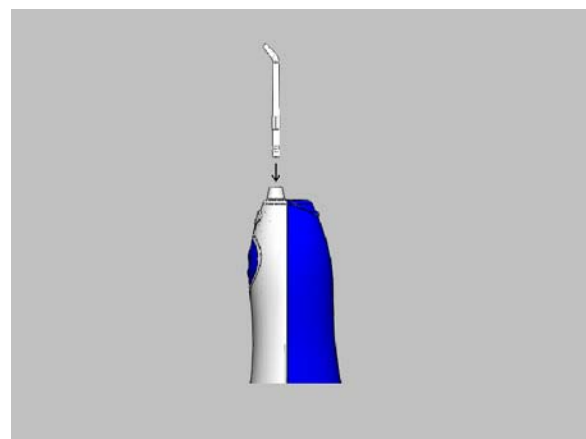**( D )**

## HOW TO USE

5. Place the Jet Tip into your mouth. Aim the tip toward the teeth, close your lips to prevent splashing while still allow water to flow from mouth into the sink, then turn the unit on ( **E** ).


**Note:** In order to prevent the water from splashing everywhere, make sure the unit is turned off before removing the Jet Tip from your mouth during use.


( **E** )

**Note:** To get the best water flowing path, please keep the unit upside up ( **F** ).

**Note:** The water stream orientation can be adjusted by turning the Jet Tip or Tip Rotator at the top of handle ( **G** ).

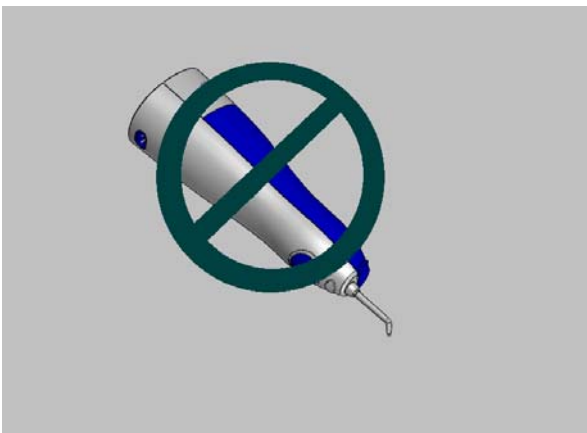
( **F** )

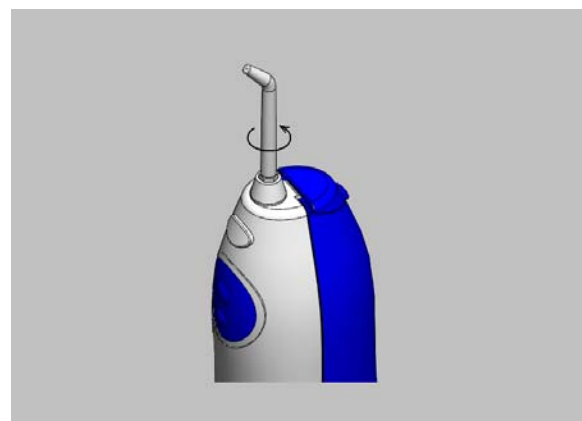
( **G** )

## HOW TO USE

**6.** Move the Jet Tip to any area between the teeth until the teeth have been cleaned ( **H, I** ).

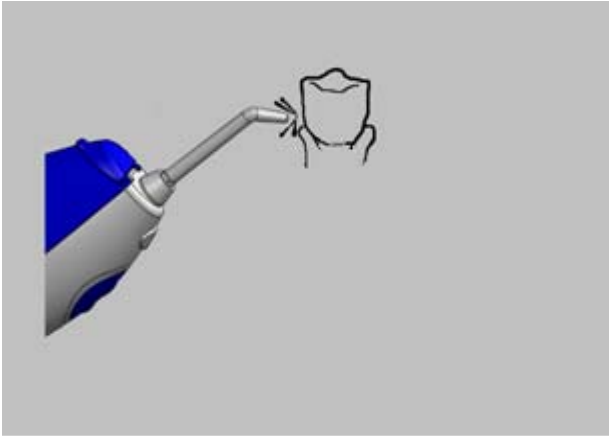
( H )

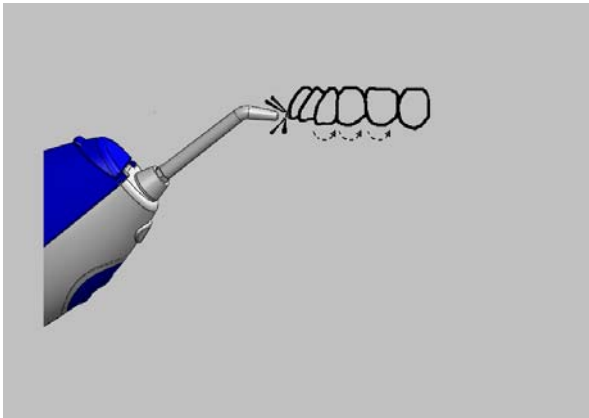
( I )

**Note:** If you want to change the Jet Tip, press the Jet Tip Ejector and the Jet Tip can be removed from the Handle ( **K** ). Refer to Step 4 for the Jet Tip insertion.

**7. When finished.** Turn the unit off. Empty any liquid left in the Reservoir ( **J** ).

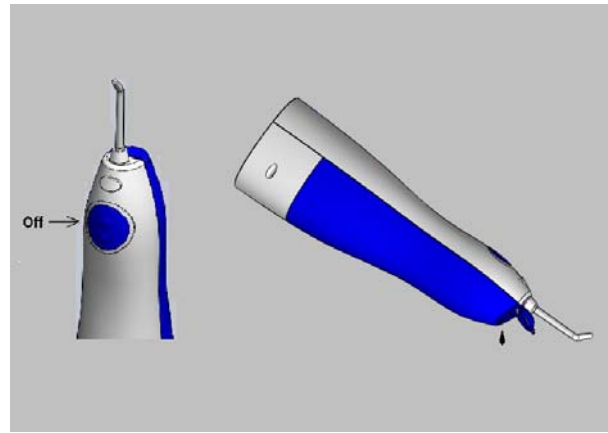( J )

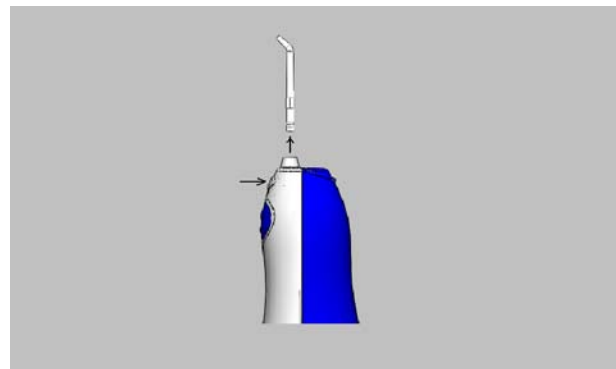( K )

## CARE AND CLEANING

### Cleaning

1. Be sure the unit is unplugged.
2. Wipe the appliance's exterior surface with a soft damp cloth.
3. Do not use abrasive cleaners. **Never immerse the appliance in water.**
4. The reservoir is not removable. If you want to clean the inside reservoir, open the reservoir lid and let the water flow in and out to clean the reservoir. **Caution: Do not soak or immerse the unit in water.**
5. Mineral deposits left by hard water or any other unrecommended solution can clog your machine. Use the clean water only or the other appropriate solution recommended by your dental professional.  
**Note: Use of some antibacterial solutions may cause the machine malfunction or shorten the life of the appliance.**

### Care

This product contains no user serviceable parts. Refer service to the authorized service center or by Candeon Technologies Limited.

### Battery Disposal

This product contains a Rechargeable Nickel Cadmium battery. The use of rechargeable batteries can reduce the number of batteries entering the waste stream but it still contains heavy metals. In order to help saving our natural resources & protection of our human health, please recycle the battery at your nearest recycle center when the unit is run out of its life. **Do not dispose battery as household waste.**

## TROUBLESHOOTING

Problem	Possible Cause	Solution
Unit does not work.	Battery out of power.	Turn Unit off. Connect unit to charger. Charge the Unit for 16 hours. The LED will light indicating the connection of unit to charger is good.
Water leaks between Jet Tip and Tip Holder.	Jet Tip not in the position.	Insert Jet Tip to the position.
No water flow.	Upside down of the unit.	Keep unit upside up.
	Run out of water inside the Reservoir.	Fill water into Reservoir.
	The flow path is clogged.	Change to another Jet Tip.
		Fill Reservoir with the clean water, add a spoon of vinegar (10ml) into the Reservoir, run Unit to empty water, fill Reservoir with water and run unit again.

## ONE-YEAR LIMITED WARRANTY

### **What does it cover?**

Any defect in material or workmanship. However, Candeon's liability will not exceed the purchase price of the product.

### **For how long?**

One year after the date of original purchase.

### **What will we do to help you?**

Provide you with a reasonably similar replacement product that is either new or factory refurbished.

### **How do you get service?**

Save your receipt as proof of the date of sale.

Contact the nearest distributor for assistance.

Check our on-line service site at [www.candeon.cn](http://www.candeon.cn) for help.

### **What does your warranty not cover?**

Damage from misuse, abuse or neglect.

Damage from commercial use.

Products that have been modified in any way.

Shipping and handling costs associated with the replacement of the unit.

Consequential or incidental damages.

**CANDEON**

Look and Feel Your Best

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**China Office**

1309 Jin Tian Building  
He Ping Road Shen Zhen, China  
Tel: 0755-2556 8355

**US Office**

10640 Merriman Rd  
Cupertino, CA 95014  
USA

**CANDEON™ and 凯电™ are trademarks of Candeon Technologies Limited.**

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